

IC factor (carb ratio)

Contribution to the discussion among DIY loopers

The author assumes no liability V.3.1 Jun24



1. What is the IC factor?

1.1 To determine a user bolus

1.2 To calculate carb absorption

1.3 IC can vary (see also 7.)

2. Rough estimates for your IC

3. „Experimental“ determinations at meal times

3.1 Basal rate first

3.2 Determinations for each meal time

3.3 Using the IC to set a meal bolus

3.4 What when your bolus wears out?

4. Determination of IC via other factors

5. Avoiding high glucose

6. IC tuning in Hybrid Closed Loop

7. Sensitivity-adjusted IC (carb ratios)

7.1 Circadian pattern

7.2 Temp. adaptations to sensitivity

7.3 Automatic adaptations using dynamic carb ratio?

8. Limited role of IC in Full Closed Loop

1. What is the IC factor?

The IC factor describes how much insulin you need (to return to your starting-bg level) for carbs that you consume.

For example, $IC = 8 \text{ g} / U$ means, for 8 g carbs, 1 unit of insulin is required

○

1.1 Usage to determine a user bolus in hybrid closed looping

The IC factor is a key parameter for your determination of a bolus via the bolus calculator (also called bolus wizard) in hybrid closed looping.

Please refer to the related instructions, like for AAPS here

<https://androidaps.readthedocs.io/de/latest/Getting-Started/Screenshots.html#bolus-wizard>

Observe that you should bolus only for as many grams of carbs as will be absorbed (digested) in the time window when your given meal bolus has really strong activity. So, either enter only up to 60, maybe 75 grams (depending on how slow your insulin is); - or you can also enter more carbs, but apply a %age to be bolussed for. See also section 3.3

Your selected insulin, the relative timing of your given bolus and meal start, and the meal composition define whether you will see one long-stretched bg hump, or (preferably, after “sharpening” your IC) an S curve, with bg going low right after the activity max from your insulin bolus, and from there, when your bolus loses power and first SMBs are needed, rising less high,.

More see section 6., and also “Meal Mgt. Basics.pdf in:

<https://github.com/bernie4375/HCL-Meal-Mgt.-ISF-and-IC-settings>.

For later absorbed carbs, the loop will take care (see sections 3.4 and 5.), without the need for you to provide extra “extended” boli.

<https://androidaps.readthedocs.io/de/latest/Usage/Extended-Carbs.html#extended-bolus-and-why-they-won-t-work-in-closed-loop-environment>

1.2 Usage to calculate carb absorption

1.2.1 Basic carb absorption calculation

(From: <https://androidaps.readthedocs.io/en/latest/Usage/COB-calculation.html>): When carbs are entered, AAPS adds them to the current carbs on board (COB). AAPS then absorbs (removes) carbs based on observed deviations to BG values according to the formula:

$$\text{absorbed_carbs} = \text{bg deviation} * \text{ic} / \text{isf}$$

For any 5 minute segment, the loop can exactly calculate the amount of insulin used up in that time segment (delta iob; e.g. *minus 0-7 U*)

For this calculation, the loop uses the kinetic insulin model, as in AAPS can be seen via the blue curve in the INS tab. See also “Insulin..DIA..pdf” in: <https://github.com/bernie4375/HCL-Meal-Mgt.-ISF-and-IC-settings>

The observed change in bg (bg delta) divided by ISF is the amount of insulin used up for the observed bg correction (e.g. *at an ISF of 40 mg/dl/U, when we see bg lowered by 20 mg/dl: 20 mg/dl divided by 40 mg/dl/U = 0.5 U*)

The remainder of **iob consumed** ($0.7 - 0.5 = 0.2 \text{ U}$ in our example) **multiplied with the IC** is the amount of carbs absorbed (e.g. *for someone with an IC = 10 g/U: 0.2 U times 10 g/U = 2 g*).

In our example (*0,7 U delta iob in 5 minutes, bg falling by 20 mg/dl and at assumed IC and ISF*), *2 g of carbs were absorbed in 5 minutes*. (This would be 24 g/hour, a very reasonable value).

1.2.2 Dynamic carb absorption

oref(1) loops do not directly use IC (and your carb inputs) to determine the automatic insulin delivery.

1.2.2.1 Carb decay (calculated, looking back)

As already explained in section 1.2.1, the loop will always **look back** at the most recent 5 minute development and use the IC (but also the ISF) **to calculate carb decay**, while assigning otherwise unlogical effects to temp. changes of your insulin sensitivity:

The algo knows exactly how much insulin is consumed in each 5 minute segment. And it sees the bg delta. From there it can simply calculate, using the ISF, how much of the consumed insulin went for bg correction. And how much, then, must have been consumed for carb decay (using the IC).

There will be instances with a not-plausible carb result (the plausibility corridor is defined by max carb absorption, usually 30g/h, and by min5m Carbimpact) , the difference is interpreted as temp. sensitivity change (->Autosens). ((It could also be ascribed to a lousy CGM, which principally must be avoided, maybe just by picking the right smoothing option)).

1.2.2.2 Expected carb decay (estimated, looking forward)

With the UAM setting ...

- even in lack of any carb inputs (I do not make any)
- regardless of what may be left of your initial meal input (g of carbs, minus absorbed carbs according to [section 1.2.1](#))
- and regardless also whether you gave inputs on future carbs (“eCarbs”)

..the algo always assumes the carb absorption that “pro forma” (in the aforementioned calculation) resulted in the past 5 minutes, **is likely to continue, then could slowly fade out over the next hour or so ...**

Background: 1) It is not possible that our digestion all of a sudden stops; it either is in a rise phase still, or at a constant “burn” level, or at a – likely a S shape – decline.
2) Your given e-Carb input might suggest less “fading out”, and this will go into one of the predictions. However, knowing how uncertain user inputs about grams of carbs coming to be absorbed in later hours after meal start are, the loop will put higher emphasis on the UAM calculation driven prediction (and, in UAM FCL, fully rely on it).

Note that, **always 5 minutes later, it can** see how it really went, and **adjust again**. This “touch and go” will very often be better than what you could attempt to tell your loop about exact grams, and when they actually will be digested.

For that reason, at least in AAPS (the author does not know the other algos in detail), also people who do make detailed carb and eCarb inputs, do NOT have their loop really run much on that inputted carb data**).

Just one of 4 predictions makes use of the carb info. And also in this case, of making carb inputs, **"carb deviations" seen by "UAM" constitute another prediction** (that, to the loop, is less doubtful than your everyday inputs).

******)If UAM is not switched on, if some of your settings depend on $cob > 0$, or if you want to use Autotune (not recommended by me) there can be merits of precise carb inputs

In conclusion, my advice for users of systems with dynamic carb absorption on board (OpenAPS, AAPS, iAPS,) would be, to put less every-day effort into getting carb estimates right, but (periodically, and especially in the beginning) more effort into determining your factors (especially "the other one", ISF), and into understanding how your algo determines SMB sizes etc.

Note 1: Tuning (factors, and dynamic factors) often is done completely wrong, if not paying attention to cut-offs (like max "allowed" basal minutes per SMB)

Note 2: Dynamic carb absorption works in principle with your profile IC. Whether there is additional merit to a "dynamic carb ratio", is discussed further below.

1.3 IC values will vary

Unfortunately, your carb ratio (IC) will not be just one fixed number you could always count on (calculate with).

1.3.1 IC varies for each meal time

The IC factor should be determined separately for each meal (breakfast, lunch, dinner). Check whether you see a typical „**circadian**“ sensitivity pattern corresponding to your bio-rhythm..(As a high basal rate and a low IC factor are signs of lower insulin sensitivity (=of elevated insulin need), the pattern must be like a mirror-image. More see in section 7.1, Circadian Pattern, below)

140 1.3.2 IC can also “situationally” (temporarily) vary

141 For instance when hormones play into it, but also for many other reasons, your
142 sensitivity to insulin may temporarily change.

143 Autosens might alert you to this, and (depending on your settings) make automatic
144 profile adjustments.

145 In some instances (like planning exercise) you would know beforehand, and should
146 manually set a timed %profile switch.

147 A % profile switch adjusts IC (as well as ISF, and also profile basal) to the observed
148 (by Autosens) or expected (by you) insulin sensitivity, see e.g

149 <https://androidaps.readthedocs.io/de/latest/Usage/Profiles.html#percentage>

150 More see [section 7.2](#)

151

152 1.3.3 Dynamic carb ratio?

153 The (controversial) hypothesis behind suggesting a dynamic IC (carb ratio) is that a
154 different IC should be applied depending on bg level, and on TDD. More see section
155 7.3

156 Having a range of ICs, automatically adjusting to likely adapted insulin needs, can be
157 good, at least in the sense of eventually, better late than never, adjusting both boli
158 and carb absorption, for some (delayed) improvement.

159 Meals should never start at super high bg, though. Dynamic carb ratio would, then,
160 never really be used for bolussing. And even to the extent there would be an
161 elevated starting glucose, the ISF would take care of the “correction” part (see also
162 bolus Calculators)

163

164

2. Rough estimate for your IC

2.1 Autotune (not recommended by the author)

Autotune gives one „average“ IC. You could make the effort to manually differentiate according to your established 24 hour pattern. Still, „reliability“ of the IC resulting from Autotune is seen controversial - to a large part probably because you don't always enter complete and 100% correct data. Inaccurate/inconsistent input => less useful output!

2.2 IC estimate based on (TDD minus profile basal)

You can get your **daily average IC** if you (1) count up the g carbs in 24 hrs and (2) divide it by the amount of 24 hr bolus insulin. Problem with the latter: Because your loop modulates basal rates all the time, you must first look up the **TDD (total daily insulin given)** (in AAPS you see that at the bottom of the /ACT/ screen, (next to HOME, or in statistics). Then subtract the "real" 24hr basal need as in your profile, from the TDD):

IC (g/U) = C (daily g carb) / (TDD - 24h Basal as in Profile).

Adult example: TDD = 37U; Profile Basal = 16U; daily carbs 200g

$\Rightarrow IC = 200 \text{ g} / (37U - 16U) = 200 / 21 \text{ g/U} \sim 9,5 \text{ g/U}$

Try to eliminate days with extreme sports, unusual stress, or infection from that evaluation. Later you will modify insulin delivery for such scenarios via profile switch = „tuning“ your IC according to the temporary changed typical insulin requirement. Therefore, avoid "averaging" such effects into your factor determination upfront already

2.3 IC estimate based on TDD

For a very rough first estimate for your IC (and other key profile factors), you can – as some commercial systems do for easing into looping – just start from TDD (as your TDD roughly describes your sensitivity to insulin, at your average diet and activity level).

193 The following is a suggestion copied from
194 <https://www.wcu.edu/WebFiles/PDFs/CalculatingInsulin.pdf>

Then for basal bolus calculate what percentage you want. Typically 40% basal and 60% bolus.

Ex) 40% of 40 units = 16 u basal & 60% of 40 units = 24u bolus total then divide by 3= 8units per meal (for 3 meals per day)

Calculating Insulin Sensitivity Factor (AKA Correction Factor) ISF

1500 divided by Total Daily Dose of insulin (TDD) if patient uses rapid acting insulin

OR 1800 divided by TDD if patient uses regular insulin

Ex) TDD = 40 units so $1500/40 = 37.5 \text{ mg/dl/U} = \text{ISF}$; divide by 18 for mmol/l/U

If current premeal BG is 160 and the target BG is 90 you would take the current BG subtract the target BG then multiply by the correction factor.

195 Ex) $(160-90)/37.5 = 1.9 \text{ units}$

Carb to Insulin Ratio IC

This is the number of grams of carbohydrates that is covered by 1 unit of insulin.

How to calculate: 500 divided by TDD $\Rightarrow \text{IC} = 12.5 \text{ g/U} \Rightarrow 8 \text{ U avg need for } 100 \text{ g carb meal}$

Ex) $500/40 = 12.5 \text{ grams per unit}$ (I:C ratio is 1:12.5)

So if 90 gram meal then you would divide 90 by 12.5 = 7.2 units

196 If target BG is above range for 2-3 days then decrease C:I ratio by 10-20%, if target BG is below range for 2-3 days then increase C:I ratio by 10-20%.

197

198 3. “Experimental” determination of IC at meal times

199

200 Attempts to define the IC factor via TDD (section 2.) can only yield a rough estimate,
201 also because the influence of fats and proteins is often omitted, or done wrong, or
202 applied inconsistently.

203 This is discussed in more detail in: “Meal-Management Basics.pdf”, see :

204 <https://github.com/bernie4375/HCL-Meal-Mgt.-ISF-and-IC-settings>

205 To determine a meal bolus in Hybrid Closed Loop, best use an IC-factor which you
206 determined for the relevant time-of-day window, as follows.

3.1 Basal rate first

Verify that you have a correct basal rate before determining factors!

This is important because you want to watch the effect of your meal bolus for at least 3 hours, and a basal rate error in that time window would bias your observations. Additionally, your observations would be burdened by any wrong basal values for a couple of preceding hours.

As a consequence, you would lose some of the principal power of your loop.

Helpful links regarding basal rate

- A basal rate helper, drawing data from your Nightscout site (by Peter van Rijt):
<https://www.nightscoutsuggestions.com/index.php>
- a good article in support of that process:
<https://www.mysugr.com/en/blog/basal-rate-testing/>
- a process description by a commercial service, including a tool :
<https://www.mevita.de/.../online-zugang-erstellung-einer.../>

Although I advise to determine your true, likely **circadian**, basal rate:

It can be „good enough“, and advantageous especially for small kids, to simplify things and assume a (in case of doubt, low!) **flat basal rate**, together with a well-contoured ISF pattern („with enough bite where it counts“). Discussions on this topic can be found in looper sites using the search term “flat basal rate”. Interpret enthusiastic reports on flat basal with caution, though. Positive observed effects are likely coming from significantly increased attention to ISF and IC... and this brings us back to our topic...

3.2 Determination of your IC for each major mealtime

Make sure your basal rate is right. On a day without preceding major activity, stress, infection, and a relative steady glucose in the normal range (and cob=0) before meal start:

Shut closed loop off. (Open loop, with just the profile basal rate running).

Eat a well defined **smaller meal** (20 .. 45 g of preferably „rapid“ carbs; not much fat and protein, please) and use your suspected (see section 2) IC to determine the amount of insulin for this meal.

With Closed Loop off, but profile basal running (i.e. in Open Loop), watch for 3 hours.

(This assumes you use Lyumjev or Fiasp, at least in a 50% mix. The author never tried any slower insulin, which would be inferior for looping, and require longer time periods to observe desired changes, in tests as well as in everyday life!)

- If your glucose levels out about where you started, the IC can be used.
- If curve goes too low (eat some carbs and) try next day again, with a higher IC value.
- If curve remains too high, the IC was too weak, and needs to be lowered.

Using your ISF, you can calculate how many units of insulin less, or more, your IC rather should have provided to actually come to bg target. Units of bolussed insulin / desired insulin would be an estimate for the factor to adjust your IC. (This is for open loop testing. In closed loop you would have to factor SMBs, + TBRs deviating from 100%, in).

For testing carb ratios see also Katie di Simone (†) at:

<http://seemycgm.com/2017/10/29/fine-tuning-settings/> and Dr. Saleh Adi from Tidepool at: <https://www.youtube.com/watch?v=McxO3oOkzc4>

3.3 Using the IC to set a meal bolus

See instructions coming with the bolus calculator (or bolus wizard) that comes with your loop or pump. For AAPS, see here:

<https://androidaps.readthedocs.io/de/latest/Getting-Started/Screenshots.html#bolus-wizard>

People who eat **very carb-rich diets** must give some consideration to the fact that the capacity of their body, how much carbs it can absorb per hour, is limited (in adults often to 30g/hour; see chapter 6.).

Also, there is a marked difference between insulins. Check (e.g. via the pink curve in the AAPS insulin tab) in which time window your bolus loses most of its activity.

For instance, for Lyumjev, at 120 minutes after injection already 75% of activity is used up.

So **only a portion of the meal might be servable via a meal bolus** given at the beginning.

The good news for high carbers on AAPS or otheroref(1) systems: You can determine your everyday meal bolus by just dividing the g carbs that are absorbable while your bolus goes strong (for Lyumjev: 60g) by your IC, and “always” bolus that.

It is a frequently seen mistake that higher carb amounts (the entire meal) are entered (to 100%) into the bolus calculator. Subsequent tuning (to avoid hypos) makes people soften up (elevate) their IC, and their loop “lacks bite” as a consequence. ((You then experience higher bg, which could be counter-balanced with - at high bg - strengthened factors (dynamicISF/dynamic carb ratio). That way many come to a solution that may often look good enough. However you easy end up creating a maze of little errors and counter-balances to them, which makes your loop unstable! - The author advises to do a solid groundwork first, before resorting to extra “tweeks”))

PS: If you see merit in announcing more carbs right upfront, then make use of the % -to-be-bolussed-for button, and keep an aggressive IC intact.

3.4 What when your bolus wears out?

The IC- driven bolus should get you through the first 2-3 hours of a meal. (More see section 6 on tuning)

Challenges can arise after the 2nd hour when meal bolusses „wear out“ and Fat/Protein contributes.´ <https://androidaps.readthedocs.io/de/latest/Usage/Extended-Carbs.html#extended-carbs-ecarbs>

Regarding late meal phase see also Meal Management in <https://github.com/bernie4375/HCL-Meal-Mgt.-ISF-and-IC-settings> or also this study: <https://www.ncbi.nlm.nih.gov/pmc/articles/PMC4454102/>

3.5 Profile helper for kids

To help “construct” a circadian profile for kids, there is a “profile helper” in AAPS: <https://androidaps.readthedocs.io/en/latest/Configuration/profilehelper.html#profile-for-kids-up-to-18-years>

4. Determination of IC ratio via other factors (CRR, CSF, ISF)

In the preceding section we learned that determining suitable IC factors is not easy.

Size and composition of a meal, the selected insulin parameters, and roundabout 40 other factors occasionally present, can introduce variations in our attempt to define our IC.

In the following, therefore another factor is presented which is much easier to determine, and can serve to either determine the IC via an alternative route, or to provide a plausibility control.

The **carb rise ratio** (by some also called CSF, carb sensitivity factor) describes by how many mg/dl our glucose rises per gram of absorbed carbohydrate.

CRR is actually quite easy to determine (easier than the IC, which requires 3 hours with underlying stable glucose a correct basal rate, and rapidly absorbed carbs):

(In a relatively stable, normo-glycemic phase) take a sweet drink with a known carb content, and watch (in Open Loop, not giving any insulin except profile basal rate) by how many (mg/dl) glucose rises, until reaching a plateau (in about 1 hour).

Example: After taking 20g of carbs, bg rises from 90 to 190 mg/dl.

$$CRR = (190-90)\text{mg/dl} / 20\text{g} = 5 \text{ mg/dl} / \text{g}$$

This parameter is of great value to check your ISF and IC values for plausibility, because:

$$\text{ISF (mg/dl drop per U)} / \text{CRR (mg/dl rise per g carb)} = \text{IC (g carb / U)}$$

or $CRR = ISF / IC$ (Likewise of course in the mmol world)

I can now continue the „experiment“ (example above):

By treating the plateau of 190 mg/dl with an amount of insulin that suits my insulin sensitivity, e.g. 2 Units. Then I observe for about three hours (ClosedLoop off, profile basal running), until a new lower plateau is reached., For instance, a new plateau might build at 110 mg/dl, in which case 2 U of insulin brought me down by 80 mg/dl. My ISF would calculate to $80 / 2 = 40 \text{ mg/dl} / \text{U}$.

$$\text{The IC would follow as: } ISF / CRR = 40 \text{ (mg/dl)/U} / 5 \text{ (mg/dl)/g} = 8 \text{ g/U} = IC$$

The CRR is relatively stable (does not vary by much over 24 hours): Both, ISF and IC, do vary (in similar ways, related to your insulin sensitivity), but when dividing on by the other, deviations owed to sensitivity cancel out.

5. Avoiding high glucose values

Even if your IC value is correct (brings you back to target about 2-3 hours after a meal), you might be disappointed by **intermediate high glucose peaks**.

Resist the temptation of extra bolussing when seeing a high peak! Rage bolussing comes with significant danger to run into a hypoglycemia a bit later.

Rather, try to picture for yourself the course of carb absorption on one hand, and of insulin activity developing, on the other.

See also “Meal Management Basics.pdf” in: (<https://github.com/bernie4375/HCL-Meal-Mgt.-ISF-and-IC-settings>).

Regarding **insulin activity**, AAPS displays that as a thin yellow line in your main screen if you select it (top right in your glucose screen press the little dart, activate Basal and Activity).

Regarding **carb absorption** it is important, that it starts before any insulin activity (hence always a rising glucose, initially), and then it can run rather steadily – in most adults at about 30 g/h as Dana Lewis has observed - Fat and fibre have an additional effect of stretching absorption out.

There are several strategies to minimize glucose peaks in the first hour or two after a meal:

1. Pre-bolussing a couple of minutes before the meal starts However, this can be dangerous if timing (when eating “must, latest” start) is not strictly adhered to
This can be problematic in restaurants, for instance, where you may need to “bridge” a delay being served ...by going for any carbs in immediate reach.
2. Give only a small part of the meal bolus before the meal begins. This enhances safety but increases complexity.
3. Orient your loop, already an hour +/- 30 minutes before any meal, towards a lower glucose goal. This strategy has 3 nice benefits: (1) It lowers the starting glucose, so the peak from the meal will be accordingly lower (2) Moreover, you get some positive job at meal start, further supporting a milder rise. (3) This move is very time-un-critical, and can even be automated for some of your meal times. ((Even if you skip a meal, nothing bad happens, other than that you need a snack in case you want to start exercise, rather than have a meal, when at the low range of your green glucose range))
4. Full Closed Loop can, under certain conditions, and after significant tuning effort, also provide solutions (for advanced loopers, only. See <https://github.com/bernie4375/FCL-potential-autoISF-research->)

6. Tuning IC in Hybrid Closed Loop


Following observations published by Dana Lewis *) , in the first 2 hours of nearly any meal, about 60g of carbs are absorbed, and – in HybridClosedLoop – a user bolus is there to largely take care. In this „main phase“ ,your **IC should be tuned such, that your glucose curve hits the 3 green (X)**:

Mastering the 3 phases of a meal

Pre- phase (1): Provide low starting glucose (X) and some active insulin when meal starts via „EatingSoonTT“

Main phase (2): An even stream of carbs, up to about 60g, get absorbed while the meal bolus is highly active. Tune your **IC**, so glucose stays in range (3* (X))

Late phase (3) Activity from meal bolus fades out, but carbs>60, protein and fat tend to drive glucose up now. ▲ SMBs, and well tuned **ISF**, control this.

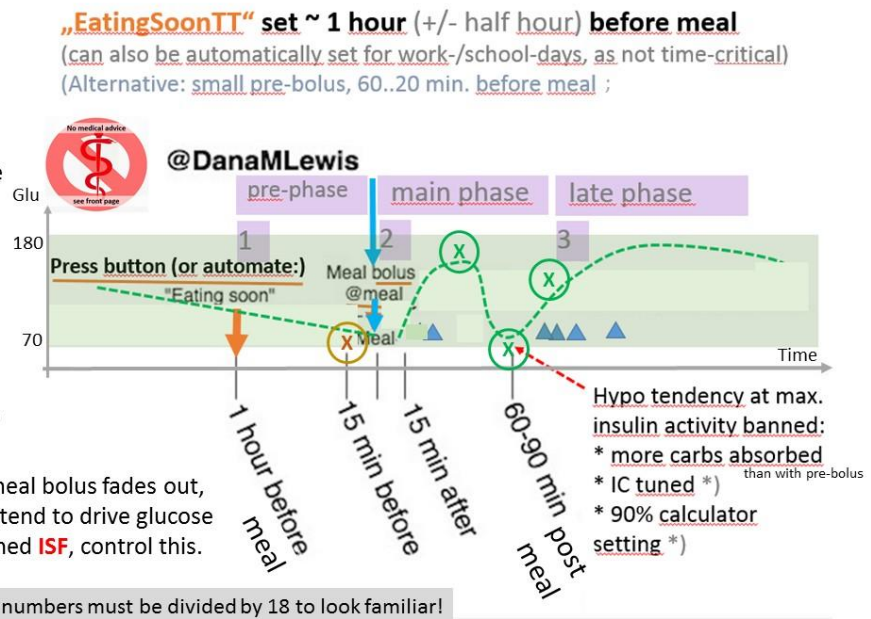

*) Chapter 8 in: <https://github.com/danamlewis/artificialpancreasbook/> -

7. Sensitivity adjusted IC (Carb Ratios)

7.1 Circadian 24 hour insulin sensitivity


Adults often have a „circadian“ pattern of insulin sensitivity „biorhythm“.

To “construct” your 24 hr circadian pattern, you can use a spreadsheet (supplied in AAPS Users FB by by Frank Duesterhoeft) if you like:

<https://docs.google.com/.../1BBOHfPFUzi4ButilakJY.../edit...> Please download the sheet and open it with Excel. All you need ist to fill the red box on Tab 1 with your values and choose a distribution Nauck, Scheiner, etc.

Your average IC, or single, in the 24 h period experimentally determined ICs, may translate into a 24 hour pattern like this:

Example

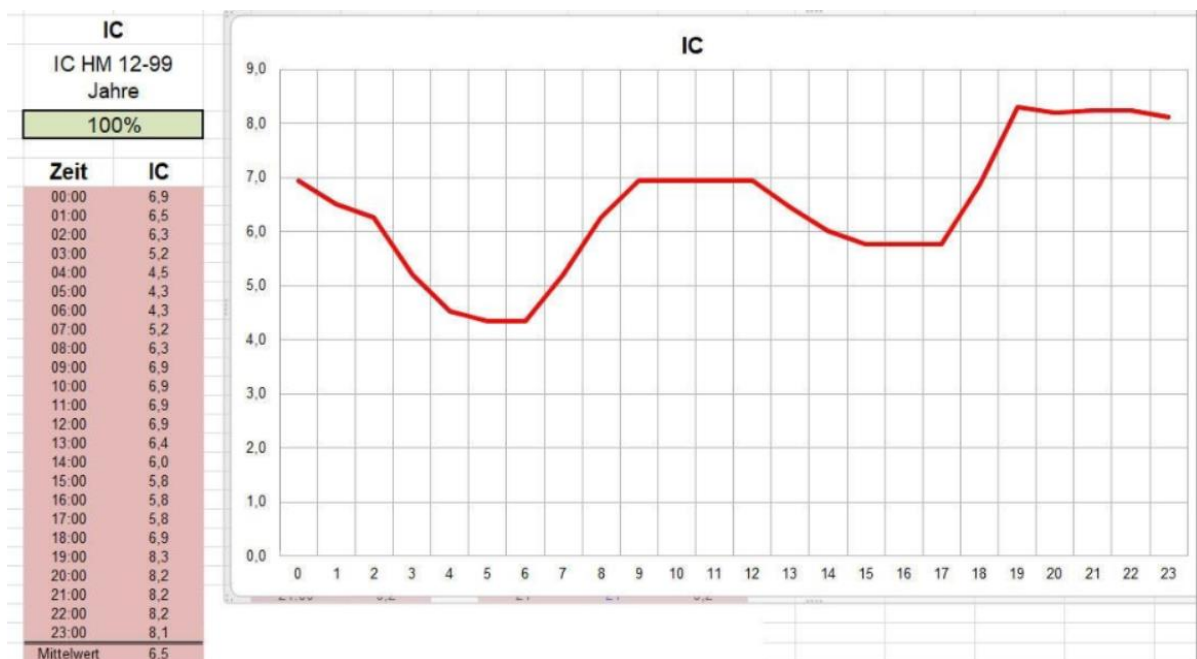

Note that lowest (most aggressive) IC is at times when your circadian basal should have highest hourly values (so, the basal curve looks inverse, mirrored on x axis).

For kids, the 24 hour sensitivity patterns can differ strongly from adults, and can also change strongly in certain growth phases. See examples given in “ISF determination..pdf” in:

<https://github.com/bernie4375/HCL-Meal-Mgt.-ISF-and-IC-settings>

To help “construct” a circadian profile for kids, there is also a “profile helper” in AAPS:
<https://androidaps.readthedocs.io/en/latest/Configuration/profilehelper.html#profile-for-kids-up-to-18-years>

7.2 Temporary effects on insulin sensitivity

Exercise , hormones, stress, (in)activity, infection (and more, see “42 factors...pdf” in:
<https://github.com/bernie4375/HCL-Meal-Mgt.-ISF-and-IC-settings>) can temporarily change sensitivity to insulin.

However, because

$$\text{InsReq} \sim (\text{Eventual bg} - \text{target}) / \text{ISF}$$

the ISF *per se* may take care of adapting to temp. sensitivity changes, while the IC would play only a minor role:

- Note, that ISF drives insulin required directly (see formula above). In contrast, IC plays only a secondary role in how eventual bg is calculated:
- The loop knows exactly how much insulin was and will be active in each 5 minute segment
- $\text{bg delta} / \text{ISF} = \text{insulin consumed for the observed bg change (delta)}$; **only the other part** of the observed delta_job is ascribed to insulin used for carb absorption:
 $\text{absorbed g carb} * \text{IC} = \text{units consumed}$
- To the extent the g carbs absorbed do not fit a corridor defined by max carb/h and by min5m carb impact, the excess is ascribed to sensitivity effects (\Rightarrow modulating IC and ISF until carb decay falls into the “allowed” corridor \Rightarrow Autosens shifts).
- The “**UAM**” part of the algo (UAM prediction) makes reasonable assumptions **how calculated past carb absorption will continue and fade out** in the next hour - and this is does **every 5 minutes** again, **rolling forward**. Very often that result is better than any user attempt at giving the loop precise info on grams of carbs, and their estimated absorption times. (Compare how your COB and UAM predictions develop after different meals!).
- g carbs inputs are not even needed at all (the loop either shows carb decay, or carb deviation), and also the IC value is not of great relevance (a lower one at high bg

would lead to slower “hypothesized” carb decay; but, as the dialed-in IC (or the cob, if any – mine is always zero, in FCL) will not really change the bg development, in essence just more will be accounted as carb deviation in each of the 5 minute segments to follow)

- => Note: g carbs entered, and also the used IC, have only little, indirect influence on how theoref loop with UAM manages glucose!)

PS: For iOS Loop, carb data are of higher relevance, and the author would be happy to add a related chapter, or good reference.

7.3 Dynamic Carb Ratio

The story started with the observation that T1Ds use stronger ISF at very **high bg** (Chris Wilson original “27700 formula”). Quickly it was noticed (through observations by Tim Street and others) that in a looping context this formula needs refinements. A personal adjustment factor, and also the observation to serve more insulin on days with **higher TDD** eventually entered the formula for dynamicISF.

- Dynamic ISF is included in AAPS Master:
<https://androidaps.readthedocs.io/de/latest/Usage/DynamicISF.html#dynamicisf-dynisf> . Note that AAPS at this time does not suggest to use dynamic carb ratio
- For a critical discussion on dynamic ISF, see e.g. “ISF determination...pdf” in:
<https://github.com/bernie4375/HCL-Meal-Mgt.-ISF-and-IC-settings>
- In iAPS, both dynamicISF and **dynamic carb ratio (IC)** are offered:
<https://discord.com/channels/1120154740857245808/1123069312295510118>

Regarding using dynamic carb ratio when determining bolus size, hybrid closed loopers should not forget the following:

- ISF is used to calculate the corrections (InsulinRequired), and the loop always applies caution in giving only 50% (then more 5 minutes later). So boosting ISF could be considered taking out some of this caution.

- However, the carb ratio IC is used mostly to determine a meal bolus = for a decision that is not every 5 minutes revisited. So, deviating much from what you had been working with in your profile seems not a good idea. (Better use the profile ICs, as were determined for each meal time)
- As it is good practice is to start meals always when bg is near target, no vast dynamic modulations of IC should happen in any case when determining a meal bolus
- Loopers should only bolus for carbs that actually get digested while their bolus is highly active. People who did that wrong in the past, had to retreat to softened IC values in their profiles. And on such basis, tightening things up now with dynamicCR can work. ((This would be a good example of two errors canceling each other out, in tendency))

As argued already in section 7.2, dynamic carb ratio might be a superfluous concept also in UAM ore loops. Often, it seems a sub-optimal “late” solution to act more aggressive at bg that already got very high. Having more aggressive ((profile, and evtl. % situational by Autosens, or manually,-adjusted)) IC and ISF already at low glucose (i.e. when bg starts getting higher) could prevent bg going high in the first place.

auto_ISF can do the best job at that.

However, there clearly are anecdotal reports from loopers benefitting using dynamicCR (and always concurrently used dynamic ISF). See <https://discord.gg/mQH9SnfeRd>

If you are interested, [Join the Dynamic ISF AAPS Discord Server!](#) for following some of the user reports on dynamic carb ratio.

To give you a glimpse into the debate there:

Chris Wilson Absent variation in meal incretin response (which governs the insulin-independent uptake of glucose from meals), the amount of insulin required to cause the uptake of a given glucose mass SHOULD vary with the base concentration of blood glucose, because the proportion of glucose disposal/uptake that is dependent on insulin scales with glucose concentration.

There are fundamental problems with the design of clamp studies that have resulted in the masking of the effects of glucose concentrations on insulin requirements. It’s an almost Schrodingerian paradox- the observation method affects the observed result.

Bjørn Ole Haugsgjerd I have also been looking at making CSF constant by scaling CR the same way as ISF. But my only reason to look at that was to hopefully improve the performance of the dynamic carb absorption model, since lowered ISF at high BGs will speed up the COB decay in ore0.

508

509 סרי נצחיה. There are definitely people that experience different insulin resistance when bg is high, among
510 them are some that are really hungry when bg is high (my kid for example) and definitely eat then... .

511

512 Too often, people praising their dynamic factors do not make the effort to investigate, which
513 part of the dynamic range was actually being used, and not falling into times of no insulin
514 required, or cut by safety features like maxSMB size (minutes of basal).

515

516 In the end you must find for yourself, which of the available methods you want to use to
517 temporarily adjust to changes in insulin sensitivity. The author finds dynamic carb ratio less
518 convincing than other methods:

- 519 • Adjusting profile% to known stages of altered sensitivity, e.g. before, during and after
520 exercise
- 521 • Using Autosens
- 522 • Defining Automations: When a certain “pattern” emerges, e.g. pointing to temp.
523 resistance from fats after a big meal, then automatically set an elevated %profile for a
524 couple of minutes (and again, as long as the condition exists).
- 525 • Using dynamicISF or autoISF

526

527 Good luck in developing your personal “good-enough” strategy for how you and your
528 loop define insulin needed for carbs. Don’t get hung up in perfectionism, enjoy times
529 with sufficiently good %TIR, as advised by your doctor.

530

531

532

533

534

535

536

8. Limited role of the IC in Full Closed Loop (UAM/no bolus/no carb inputs)

With detailed carb (amounts + absorption times) inputs, the loop has best-possible info to provide „the best expert fit“ of insulin activity and carb absorption.

It still rarely can come close to physiological values, because the time-delays inherent in our „artificial pancreas“, notably the stretched out DIA, make it difficult still, compared to a real pancreas.

So, *precise* carb inputs *could* help. However, only to the extent amounts *and* time pattern for absorption („eCarbs“) are correct ((which, every day, is a mission impossible))

Entering *imprecise* carb info could easy be inferior to not doing *any* carb inputs = to letting the *UAM mode* of oref(1) figure out further carbs that probably come to be absorbed in the next minutes, judging from the pattern of the calculated past *carb deviations* (see section 1.2.2 and <https://openaps.readthedocs.io/en/latest/docs/While%20You%20Wait%20For%20Gear/Understand-determine-basal.html#understanding-the-basic-logic-written-version>).

PS: Because that is so, also loopers *who do carb inputs* get the UAM predictions besides their other predictions, and their algo makes a judgement (every 5 minutes) as to what the best calculation might be for where glucose, underlying „real“ carb absorption, and estimated carb deviation are headed.

In any case, well-determined IC values are nice-to-have also when using the dynamic carb decay model of oref(1), also when not making carb amount entries.

Full Closed Looping

Currently, looping without carb inputs and without giving a user bolus is only possible with good results when using the oref(1) algorithm's SMB+UAM feature, as offered by OpenAPS, AndroidAPS and iAPS. ((iOS Loop, in contrast, requires fairly exact carb inputs)).

In Full Closed Loop (UAM/no bolus/no carb inputs), ISF (not IC) is the key factor for the loop to keep glucose in range.

566 The **IC factor** plays only a **minor role** there. However, the loop still uses it „in a side role“ for
567 calculating deviations = to conclude how many carbs „must have been absorbed“ in each
568 **past** 5 minute segment.

569 Note that the UAM Full Closed Loop is *not clueless* regarding how carb absorption
570 will go on, even if you did *not* give it any “extended carb” entries, and your loop
571 stubbornly sits at cob=0 all the time:

572 It will work with a **prediction** of further carb absorption building on the **carb**
573 **deviation** (=hypothesis of how much got absorbed in the past 5 minute segments),
574 and phase out more carb decay in the course of the next 1 to max 3 hours.

575 This was already discussed in section 1.2.2 on dynamic carb absorption.
576 For more detail see
577 [https://openaps.readthedocs.io/en/latest/docs/While%20You%20Wait%20For](https://openaps.readthedocs.io/en/latest/docs/While%20You%20Wait%20For%20Gear/Understand-determine-basal.html#understanding-the-basic-logic-written-version)
578 [%20Gear/Understand-determine-basal.html#understanding-the-basic-logic-](https://openaps.readthedocs.io/en/latest/docs/While%20You%20Wait%20For%20Gear/Understand-determine-basal.html#understanding-the-basic-logic-written-version)
579 [written-version](https://openaps.readthedocs.io/en/latest/docs/While%20You%20Wait%20For%20Gear/Understand-determine-basal.html#understanding-the-basic-logic-written-version) (or study your SMB tab info).

580

581 This UAM prediction about further carb absorption can be worse, but can also be
582 better than a prediction based on the user's „e-Carb“ input in Hybrid Closed Loop.

583

584 In any case, and even when having perfect knowledge about how exactly the carbs
585 fade out in the next hours, there would still be a principal problem for any loop:

- 586
- 587 • Heavy insulin „fire“ against highs will not work immediately (depending on the
588 insulin's time-to-peak), and notably it comes with a significant hypo danger
(from the „tail“ of insulin activity.)
 - 589 • A big bolus, or even a series of boli, will rarely work for several hours matching
590 the absorption of carbs (from what, how much and how fast the user ate).

591

592