

Settings for your insulin (DIA and time-to-peak)

Contribution to the discussion among DIY loopers

The author assumes no liability

B/ V.2.1 Jan'24



-
1. Setting insulin related parameters
 - 1.1 Mathematical model used
 - 1.2 Time-to-peak and DIA
 - 1.2.1 Insulin choice matters for profile ISF, IC
 - 1.2.2 Duration of insulin action
 - 1.2.3 Quantitative effects of changing DIA
 2. Other factors of potential relevance
 - 2.1 Age (of the diabetic)
 - 2.2 Dose
 - 2.3 Scatter (imprecision)
 3. Mixes of two insulins
 4. U200 insulin
-

Before doing any other tuning, make sure you have reasonably set the insulin-related parameters for your looping system.

In case you are just starting to loop and need to “household” with your time, all you need from this paper should be just two messages:

- Select your insulin in AAPS configuration, and refer to the data given in [section 1.2.2](#) regarding DIA. To set it on 7 h is a fair guess for a start, if you are uncertain.
- You probably heard that Lyumjev or Fiasp are in principle best for looping. In [section 1.2.1](#) it is explained why it is a good idea to switch at the start of your looping journey, rather than at some later point.

Changes between insulins with similar time-to-peak, like Fiasp -> Lyumjev, will be easier, and not require much of an adjustment as the case in the example given in [section 1.2.1](#).

But of course you can switch at any later time, as well. Many prefer actually to start looping with a less reactive insulin. This enhances safety in your initial months of getting to know, and tuning, your loop. (Same thought is behind the Objectives in AAPS that give you access to SMBs only after a couple of other steps).

Besides time-(minutes) to-peak activity, also the DIA (hours) selected in your profile strongly influences how the loop calculates the activity from insulin, as it unfolds in every 5-minute segment that your loop analyzes.

1. Setting insulin related parameters

1.1 Mathematical model used

Especially what should be selected as duration of insulin action (DIA) is very strongly influenced by the model used to figure out active insulin two, three, and more hours after administration. Misunderstandings about this is often a source for disputes between loopers and their treating physicians.

All insulin administrations (bigger and minor) add up to a insulin activity pattern. In the case of looping, with user boli, basal insulin, TBR modifications and SMBs given at various times, with overlapping DIAs, this can be quite complex.

In AAPS you can see insulin activity in your main screen as an extra thin yellow curve. Together with carb absorption is “explains” most of what you see in your glucose curve.

This **insulin activity pattern is an extremely important basis for each of your loop's decisions**. Having the wrong settings would give your semi-automated insulin management a permanent drift towards over- or towards under-corrections.

The loop system can still counter-regulate, but – if you burden your's with wrong DIA or time-to-peak settings in your profile – this would “use up” some of it's (limited) capacity to regulate for you.

Example: After heavy dinner, a DIA set too short “tells your loop” that active insulin is practically gone after time X. The loop takes that info for granted, and if it sees some insulin needed at that time X (and be it only for your profile basal need - as you also communicate to the loop, you need to remain stable -), then, at night-time, the loop will give you more insulin than you really need.

Therefore, before you tune your ISF differently, make sure to have a look at your DIA setting.

69

70 Please understand (and see to it, that your treating professionals understand) that models
71 can differ strongly:

- 72 • DIY looping systems use the – less common – exponential decay model.
- 73 • Medtronic uses non-linear capped curves (as in handbook to their pumps)
- 74 • Doctors / diabetes educators mostly have a rough linear model in mind
- 75 • xDrip uses a bilinear math (“with kinks”) to model insulin activity (Caution: This info
76 might be outdated)

77 All models are working „good enough“ for their (main) intended applications. But, as
78 explained above, it is worth the effort to use an exact modelling of insulin activity for a loop,
79 so it can perform optimally.

80 As pointed out already in the section 1 headline, and further explained below, the
81 mathematical model of insulin activity over time anchors on time-to-peak (minutes) and on
82 DIA (hours) in characteristic ways. This is quantitatively shown for exponential decay models
83 in [section 1.2.3-](#)

84 In AAPS, the insulin tab shows two curves: The pink one starts at 1.0 (100%) and goes down
85 to 0 (0%) when the DIA is over. It shows how much of the total activity (the capacity to lower
86 bg) is left, at any time. The blue one shows how the activity goes high, and then fades out,
87 over the DIA period (with a maximum at time-to-peak).

88

89 1.2 Time-to-peak activity and DIA for various insulins

90


91 Principally, there are „correct“ settings to each insulin type, notably regarding time-to-peak
92 activity. This is pre-programmed in the insulin choices for AAPS, for instance.

93 In the following, mostly data David Burren published in bionic wookie are cited or
94 summarized.

95

96 1.2.1 Insulin choice matters for profile ISF, IC

97 The following chart is *the inverse* of the pink curve in the AAPS insulin tab: *Not insulin still*
98 *there to be used*, but Insulin used up, going from 0% towards 100% in the 9 h DIA, for
99 Humalog with 75 minutes, and for Lyumjev with 45 minutes time-to-peak.

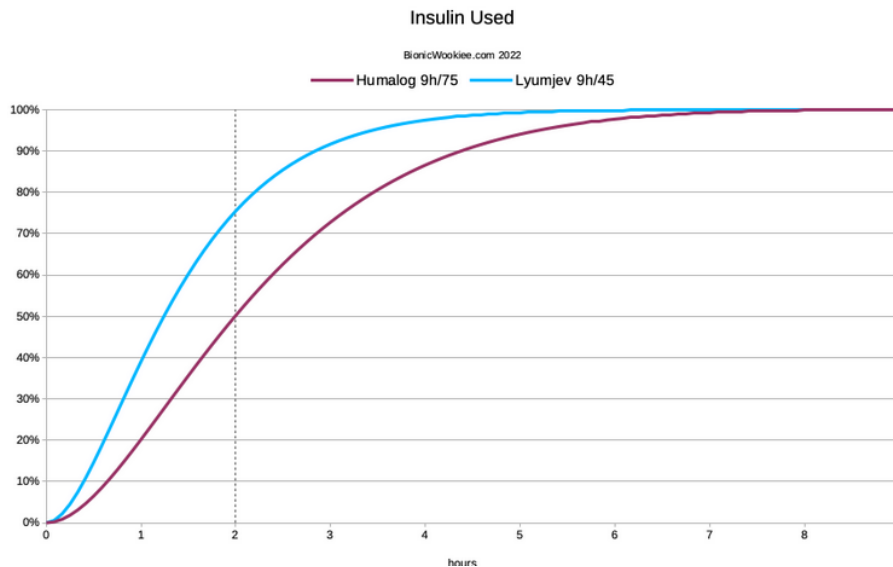

From a simplistic point of view, you can see that at the two-hour mark, more of the Lyumjev (75.5%) should have had effect than the Humalog (50.2%).

So when we're calculating *how much insulin to give for a correction*, we should tell it to give more Humalog up-front to get the same result after 2 hours.

The system will of course be tracking the IOB and forecasting the BG curves for hours into the future, so we do have some safety built in regarding the extra insulin.

When going from Humalog to using Lyumjev, this must have some consequences for the Insulin Sensitivity Factor (ISF) to use in the profile. If, for example, you had 1.8 mmol(I/U for Humalog, you should expect a "good ISF for going with Lyumjev" in the area of 2.7 mmol/I/U. According to the curves shown above (at dotted 2 hr line) a factor 75.5/50.2 applied yields the same amount of insulin for a correction.

Likewise, the Carb Ratio (IC) may deserve an adjustment when switching insulins.

- The IC could be adjusted by the same factor, for instance it might go from 7.7 g/U (Humalog) to 11.6 g/U (Lyumjev).

For a meal of 60 g, 7.8 U (=60/7.7 g/U) Humalog would have contributed 3.9 U (=50.2%*7.8 U); likewise, 5.2 U (=60g/11.6 g/U) would have contributed 3.9 U (=75.5%*5.2 U)

- For meals bigger than about 60 g you should observe that, while your insulin bolus has good activity, only a limited number of carbs can get digested (30 g/h seems the limit for most).

120 Refer to the paper on IC determination (section “Determination at meal times”) in:
121 <https://github.com/bernie4375/HCL-Meal-Mgt.-ISF-and-IC-settings/tree/FCL-w/autoISF>

122 The given example showed that switching to a “faster” insulin can have relevant
123 consequences for your key profile parameters.

124 David Burren also reports that between the two rather extreme insulin choices he tested, the
125 total amounts of insulin (TDD) did not significantly differ (- as we would expect: The same
126 amounts just gets delivered slower, even at same selected DIA, with Humalog).

127 But while the TDD has *not* changed, the instantaneous insulin levels *have*.
128 When the system is fighting post-meal “highs” the IOB will be noticeably lower with Lyumjev.

129 Although the average overall level remains similar, this might have some implications
130 for the concept of hyper-insulinaemia. This may be a subtle advantage of faster
131 insulins.

132

133 1.2.2 Duration of insulin action

134 The following focusses on the more uncertain topic of which duration of insulin action (DIA)
135 to use. It largely relies on, and quotes, results from several thorough investigations done
136 by David Burren: (<https://bionicwookiee.com>)


137 The numbers he ended up preferring are:

Insulin	Peak	Duration (DIA)
Humalog	75 minutes “Rapid-Acting Oref”	7 hours *)=
NovoRapid	75 minutes “Rapid-Acting Oref”	9 hours
Fiasp	55 minutes “Ultra-Rapid Oref”	9 hours

138 *) Later investigations <https://bionicwookiee.com/2022/04/13/revised-humalog-model-in-a-closed-loop/>
139 led to suggesting 9 hours DIA also for Humalog

140 The default constraints in AAPS have the duration limited to 7 hours, so he had to
141 make some local changes to the limits. It’s also possible if you set your “patient type”
142 to “Pregnant”, but if so you need to carefully check all the affected safety limits ([listed](#)
143 [in the AAPS documentation](#)). This may change in a future update to AndroidAPS.

For Lyumjev (45 minutes; Lyumjev Oref), there is not a big difference between a 7 and a 9 h DIA:

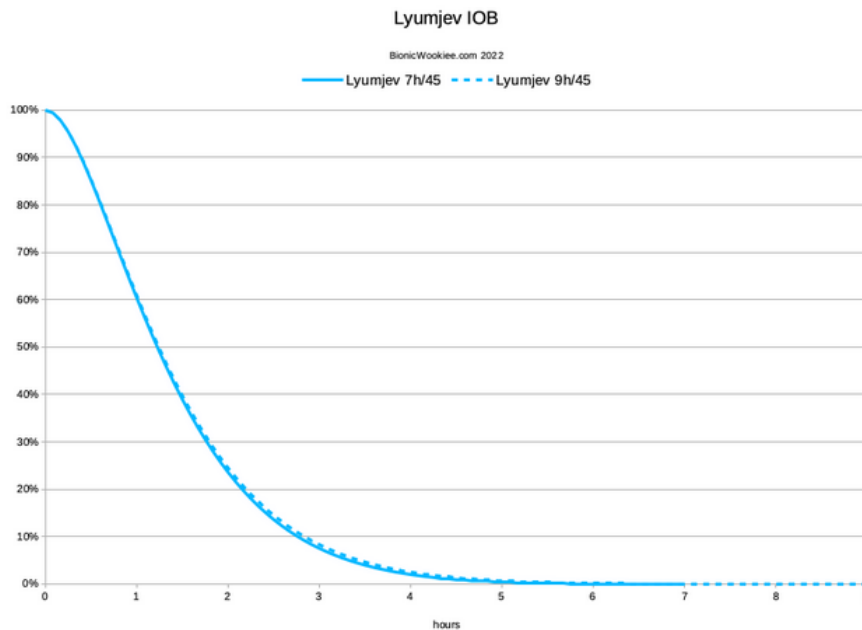

However, David Burren (<https://bionicwookiee.com/2022/04/13/revised-humalog-model-in-a-closed-loop/>) reports that, despite it's a very subtle change, he has found it can make a significant difference 5-6 hours after a meal, when ...the tails of the earlier doses do add up, and the system had been underestimating the IOB when calculating (using the shorter DIA) what was needed with new doses. With changing to a longer DIA, his average Time Below Range has reduced. (Comment: This is interesting for "fine-tuning extremists" but probably only a formal gain of little clinical significance, assuming the bg curve just swinging a bit more often, or longer, by a few mg/dl below the 70 mg/dl cut-off, that defines "below range". Judge from your own data, when/if fine tuning.)

On the DIA topic for various insulins see also: <https://www.diabettech.com/insulin/why-we-are-regularly-wrong-in-the-duration-of-insulin-action-dia-times-we-use-and-why-it-matters/>

1.2.3 Quantitative effects of changing DIA

Any given insulin dose comes with a defined total capacity for a certain bg lowering effect. How strong or weak this unfolds over a couple of hours can be mathematically modelled. In oref(1) systems, time-to-peak and DIA completely define this curve.

We can look on effects of increasing the set DIA in terms of how insulin activity would differ at any moment after administering a dose.

The next example given (chart below) does that for going from a 5 h DIA, a 6 h DIA or a 7 h DIA towards 9 h for Lyumjev

We see the peak going lower, and the tail activity higher when DIA is increased:

LYUMJEV peak @45m max effect on "tail" at ~ 3.5 h (220 minutes)

DIA 5 → 9h | **minus 5.5 %** **plus 2.5%**


DIA 5 → 6h | minus 2,7 % plus 1.3%

DIA 6 → 9h | **minus 2.8 %** **plus 1.2%**

DIA 6 → 7h | minus 1,4 % plus 0.6%

DIA 7 → 9h| - **minus 1.4 %** **plus 0.6%**

So, the „tail“ effects differ by less than 3 percent (of peak activity=100%) in the later stages of DIA:

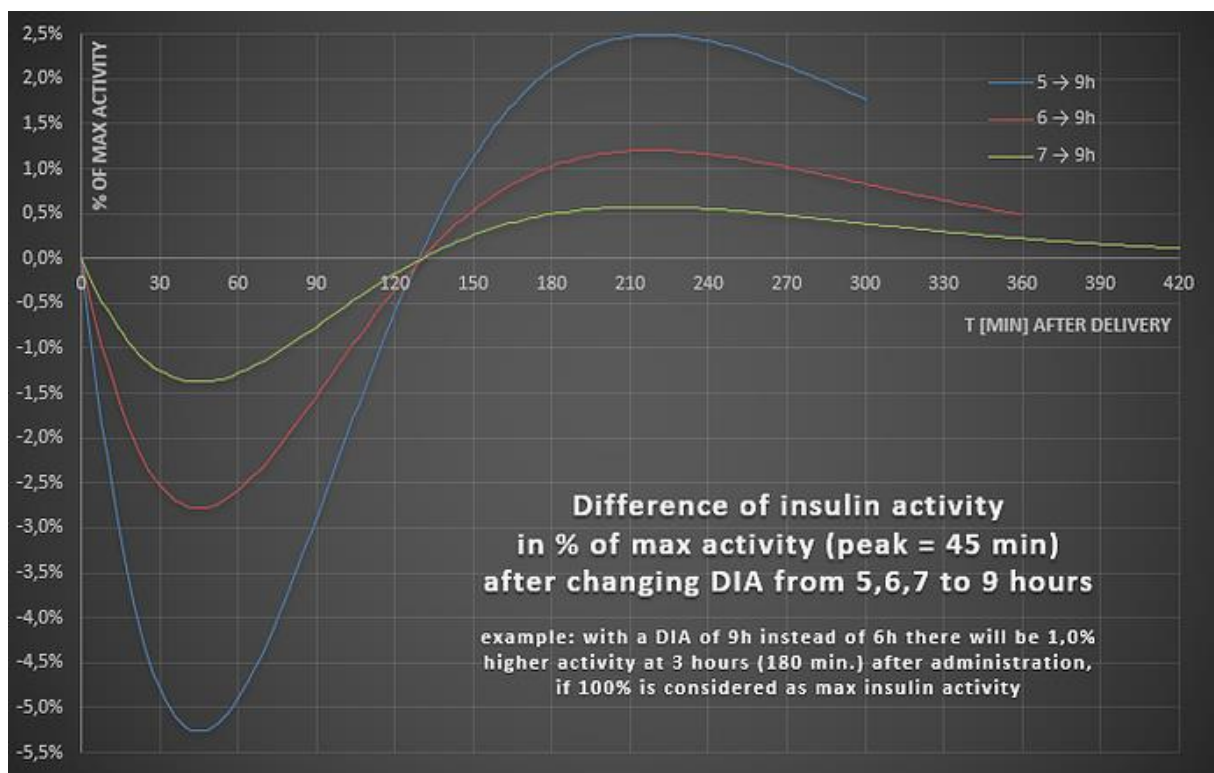

While 3 % sounds low, the significance of the problem should not be underestimated:

- For our Lyumjev case, note that the quoted 3% result is 3% of maximal activity.
Example: Activity at 180 minutes is about 0.0010 compared to 0.0080 at peak (blue curve in AAPS INS tab). 2.5% of 0.0080 would be 0.0002. BUT: 0.0012 is 20 % more than 0.0010, so **REALLY the difference in activity at 180 minutes can be up to 20%** . Still, after a bolus of 8 units (and/or SMBs that reach that iob level) for a

typical meal, the max. difference from 5 -> 9 hour DIA would roughly be, whether 1.0 U or 1.2 U are active iob left at 180 minutes. That **difference (+ 0.2 U)** should be within the loop's regulating capacity from reducing basal.

However, it becomes much bigger for users of other insulins (with longer time-to-peak):

- The delta **effects get much bigger with insulins that have a longer time-to-peak**
Some quantitative data for other insulins are as follows:

FIASP (peak=60m) min/**max** differences
DIA 5 → 9h | 6 → 9h | 7 → 9h: -10,1 / **+6,8%** | -5,6 / +3,0% | -2,9 / +1,4%

NOVORAPID (peak=75m) min/**max** differences
DIA 5 → 9h | 6 → 9h | 7 → 9h: -15,4 / **+14,1%** | -9,1 / +7,0% | -4,8 / +3,0%

*Above example applied to Novorapid **):* The effect would be up to +14.1% of max (!)
=> 2.1 U instead 1 U at 180 minutes. A **difference of + 1.2 U** results here, if DIA is set at 5, not at 9 h, so **REALLY** it could go **up to + 120% !)**

More see: szantos, de.loopercommunity.org May 2022

<https://de.loopercommunity.org/t/naechtlicher-unterzucker/10626>

**) 2,5% → + 0.2 U ergo 14.1% → +1,1 U stimmt insofern nicht ganz genau, als man beim Novorapid Case auch die Novorapid Peak-Höhe zugrunde legen müsste (die ich aber nicht greifbar habe). Wenn diese von Haus aus 20% niedriger nur kommt, hätten wir ca + 0,9U, also weiterhin etwa eine Verdoppelung ... die wir mit unserer Wahl eines längerem DIA unserem Loop sagen könnten, damit er entsprechend weniger zu-schießt ... ergo weniger Hypogefahr hinten heraus ...

2. Other factors of potential relevance


Source: [szantos](#)

The findings reported below can give you hints in which direction to look if you attempt to fine-tune your settings further, from the standard suggestion what should be suitable for your insulin ([section 1.2.2.](#)).

2.1 Age (of the diabetic)


ema.europa.eu

[novorapid-h-c-258-p46-0044-epar-assessment-report_en.pdf 3](#)

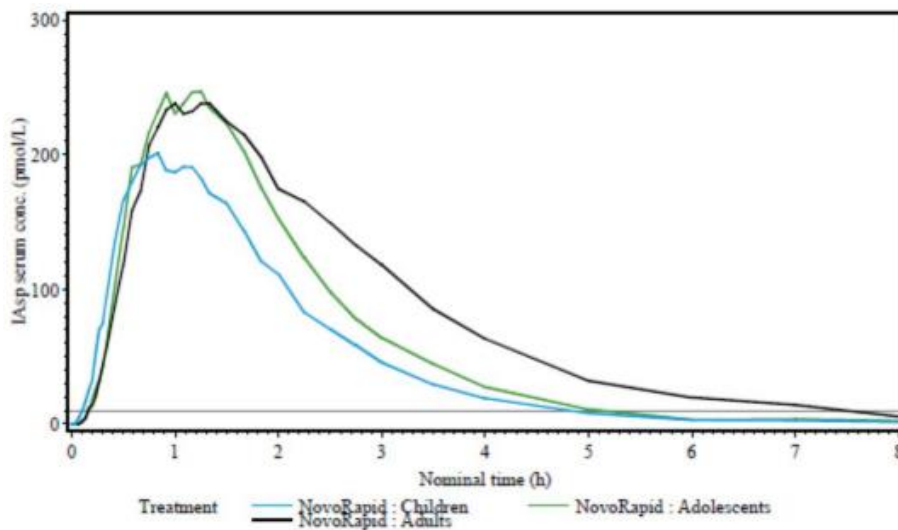

2.2 Dose

<https://journals.sagepub.com/doi/10.1177/1932296813514319>

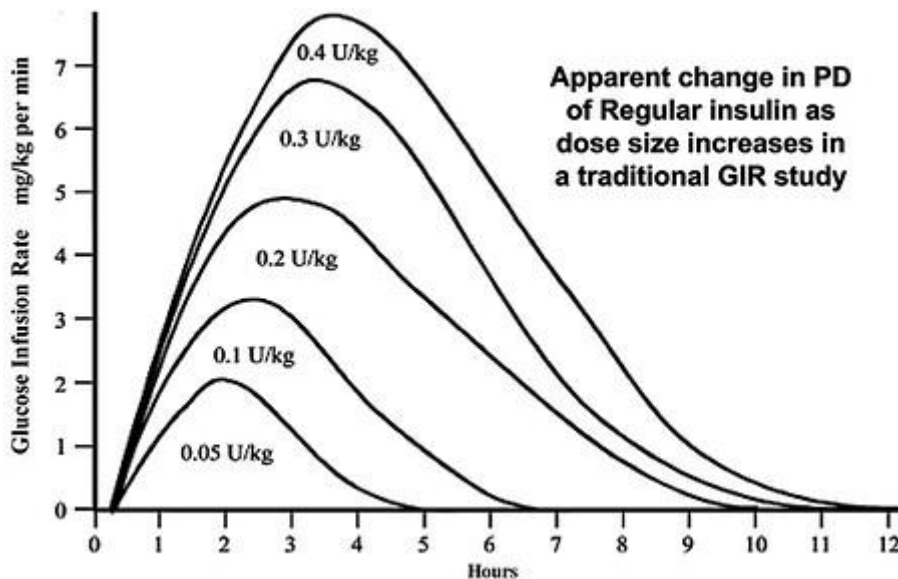

[image1000×632 78.4 KB](#)

226

227 2.3 Scatter (imprecision)

228 Individual deviations from standard suggestions could also be justified by the fact, that all
229 studies that underly the previously reported suggestions, come with very significant person-
230 to-person scatter.

231

232 All lines in the charts, as above shown from studies, are averaged data. (Some studies are
233 indicating the very significant scatter seen, as well).

234 [https://www.researchgate.net/figure/Pharmacodynamic-profiles-Insulin-action-as-expressed-](https://www.researchgate.net/figure/Pharmacodynamic-profiles-Insulin-action-as-expressed-as-GIR-required-to-maintain_fig1_41424712_2)
235 [as-GIR-required-to-maintain_fig1_41424712_2](https://www.researchgate.net/figure/Pharmacodynamic-profiles-Insulin-action-as-expressed-as-GIR-required-to-maintain_fig1_41424712_2)

236

237 3. Mixes of two insulins

238

239 The author did for some time successfully use a 50/50 mix of Fiasp and Novorapid, applying
240 the time-to-peak for Fiasp, and longest of the two DIA, as was suggested at the time, for
241 these insulins.

242

243 For a more thorough discussion see [https://bionicwookiee.com/2022/03/02/mixing-insulins-](https://bionicwookiee.com/2022/03/02/mixing-insulins-theory-and-practice/)
244 [theory-and-practice/](https://bionicwookiee.com/2022/03/02/mixing-insulins-theory-and-practice/)

245 and also: <https://bionicwookiee.com/2023/06/03/arcane-lyumjev-experiments/>

246

247

248

249 4. U200 insulins

250

251 Using up-concentrated insulins, e.g. in a U200 form, is sometimes chosen by loopers

- 252 • to reduce needed daily insulin volumes and get more time from 1 pump filling (pod)
- 253 • to reduce volume per injection for getting better tolerance regarding occlusions or
- 254 pain

255

256 There are no relevant effects on insulin parameters like DIA and time-to-peak.

257

258 However, dilution or up-concentration factors are highly relevant for setting profile factors like
259 ISF and IC, and also for some important safety settings like max iob for instance.

260

261 Refer to special discussions on that topic, e.g. here re. U200 Lyumjev

262 <https://www.diabettech.com/lyumjev/living-with-lyumjev-almost-a-year-in-review/> :

263 and also: <https://bionicwookiee.com/2023/06/03/arcane-lyumjev-experiments/>